



SURGICAL TABLE PADS

Designed for pressure management and better performance, Skytron surgical table pads feature multi-layer technology and leak-proof seams. The high-density base foam and top layer of soft, cradling memory foam provide additional levels of support and comfort.

Pads are available in a variety of thicknesses, and in bariatric sizes.



WELDED SEAMS

Unlike standard pads which allow liquids to enter between stitches and into the foam layers, Skytron pads are tightly sealed for cleanliness and durability. Ultrasonic welders are used to form a leak-proof barrier, 'welding' the seams and preventing liquids from entering the inner foam.



Welded Seam Pad



Standard Stitched Pad

MULTI-LAYER TECHNOLOGY

Skytron table pads are composed of three layers of foam, providing optimal comfort and equalized pressure regardless of patient size or procedure length. Soft, body forming, layers are reinforced with different density foams to provide lasting support.

1" viscoelastic memory foam
to equalize pressure

1/2" high density viscoelastic
foam to provide long lasting
support and elasticity

1/2" vinyl nitril closed cell
foam to strengthen the softer
layers above



AIR VENT

Air enters and leaves the pad via a covered vent, allowing the foam to expand and contract with the patient while keeping fluids out.



SURGICAL TABLE PADS

PART NUMBER IDENTIFICATION

Part numbers for Skytron pads and pad sets follow a consistent format. The numbering system always starts at the head end of the table and moves towards the foot end. In case of side-by-side pads, such as the split leg attachment, the first number will always be on the patient's left. Part numbers for pad sets will only have seven characters, while part numbers for specific pads will have eight.

PAD LOCATION



EXAMPLE

- 1 - 011 - 03 - 2 3003 table, 2" thick, pad set
- 1 - 011 - 05 - 43 3500B table, 4" thick, leg section
- 1 - 011 - 01 - 24 6000 table, 2" thick, split leg, patient's right

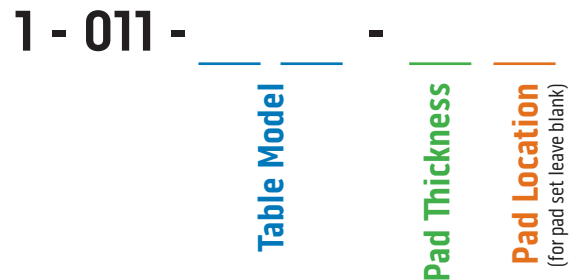


TABLE MODEL

01	1200, 5000, 6000
02	300, 1201(K), 5001, 6001(K), 6002, 6002B, 6500, 6500HD, 6600 6600B
03	3000, 3001, 3002, 3003
04	3100
05	3500, 3500B
06	3501B, 3501C, 3502, 3600 series
07	1700
08	6700 Series
09	5601E
11	6300 series
12	3501B, 3501C, 3502 and 3600 series bariatric pad set with 5" side extension
13	6700 series bariatric pad set, with 5" side extension
15	1602

The literature has been developed to provide healthcare professionals with general product information. Please note: not all products are available in all countries. Product labeling and associated claims may differ based on regulatory requirements. For product availability or any further information, please contact your local Skytron representative.



5085 Corporate Exchange Blvd. SE
Grand Rapids, MI 49512
1.800.SKYTRON (759-8766)
email: info@skytron.us
www.skytron.us

About Skytron Skytron is a privately held company specializing in capital equipment that drives efficiency in today's healthcare facilities. Our surgical tables, surgical & exam lighting, booms, video integration systems and other products help hospitals around the world work better. To learn more, visit our website at www.skytron.us.