



Skytron

Surgical Table Pads

Welded Seams

Unlike standard pads which allow liquids to enter between stitches and into the foam layers, Skytron pads are tightly sealed for cleanliness and durability. Ultrasonic welders are used to form a leak-proof barrier, 'welding' the seams and preventing liquids from entering the inner foam. Cover material offers four-way stretch to reduce damage from shear force.



Welded Seam Pad



Standard Stitched Pad


Multi-Layer Technology

Skytron table pads are composed of three layers of foam, providing optimal comfort and equalized pressure regardless of patient size or procedure length. Soft, body forming, layers are reinforced with different density foams to provide lasting support.

1" viscoelastic memory foam to equalize pressure

1/2" high density viscoelastic foam to provide long lasting support and elasticity

1/2" vinyl nitril closed cell foam to strengthen the softer layers above

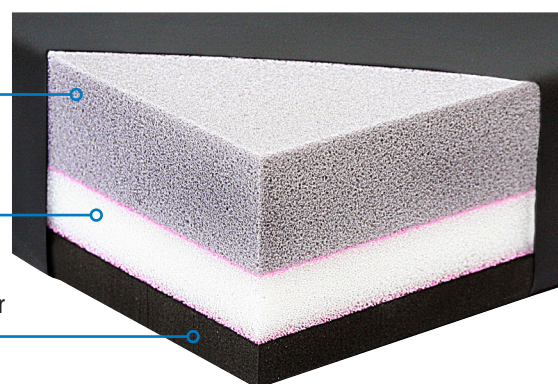

Air Vent

Air enters and leaves the pad via a covered vent, allowing the foam to expand and contract with the patient while keeping fluids out.




SKYTRON[®]
Reliability • Versatility • Cost Savings

Skytron


Surgical Table Pads

Part Number Identification

Part numbers for Skytron pads and pad sets follow a consistent format. The numbering system always starts at the head end of the table and moves towards the foot end. In case of side-by-side pads, such as the split leg attachment, the first number will always be on the patient's left. Part numbers for pad sets will only have seven characters, while part numbers for specific pads will have eight.

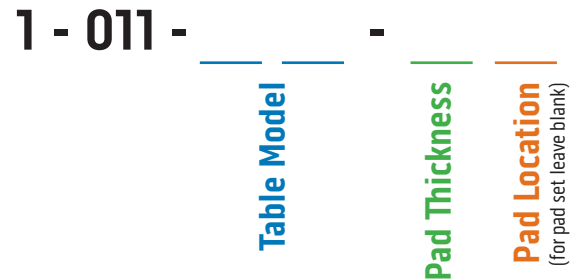

Pad Location

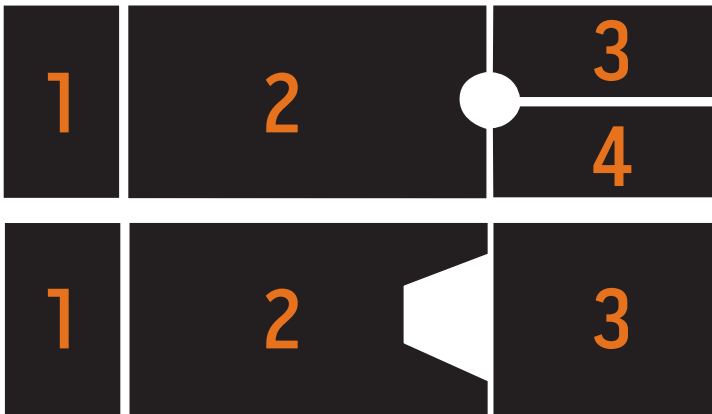

Example

- 1 - 011 - 03 - 2 3003 table, 2" thick, pad set
 1 - 011 - 05 - 43 3500B table, 4" thick, leg section
 1 - 011 - 01 - 24 6000 table, 2" thick, split leg, patient's right

Table Model

01	1200, 5000, 6000
02	300, 1201(K), 5001, 6001(K), 6002, 6002B, 6500, 6500HD, 6600 6600B
03	3000, 3001, 3002, 3003
04	3100
05	3500, 3500B
06	3501B, 3501C, 3502, 3600 series
07	1700
08	6700 Series
09	5601E
11	6300 series
12	3501B, 3501C, 3502 and 3600 series bariatric pad set with 5" side extension
13	6700 series bariatric pad set, with 5" side extension
15	1602

About Skytron: Skytron is a healthcare efficiency company specializing in capital equipment that drives performance in today's healthcare facilities. Our solutions enhance the utilization of people, facilities and capital because they are designed with the user in mind and have a low, long-term cost of ownership. We offer full-room solutions in Clinical, Infection Prevention and Clinical Business Intelligence products.

To learn more, visit our website at www.skytron.com



5085 Corporate Exchange Blvd. SE
 Grand Rapids, MI 49512
 1.800.SKYTRON (759-8766)
 email: info@skytron.com
www.skytron.com



REV5 05-19
 PN L2-010-39