




5085 Corporate Exchange Blvd. S.E.
Grand Rapids, MI 49512 • 616.656.2900

UNIVERSAL SIDE RAIL CLAMP II

P.N. 6-040-12

INSTRUCTIONS FOR USE:

The SKYTRON Universal Side Rail Clamp II can easily be attached and removed from a SKYTRON table with a U.S. standard side rail. Align the Universal Side Rail Clamp II above the side rail and lower into place. Insert the desired accessory and tighten the handle to lock into place. The SKYTRON Universal Side Rail Clamp II is designed to fit Round 5/8" (16mm) accessories.

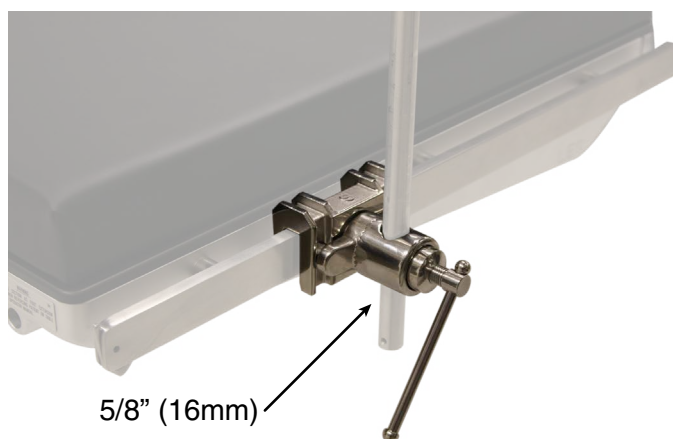

Replacement Parts:


Handle P.N. 6-040-12-1

Typical Dimensions:

Width (W) 3 1/4" (82mm)
Length w/Handle (L) 5" (127mm)

Weight:

1.95 lbs. (.88kg)

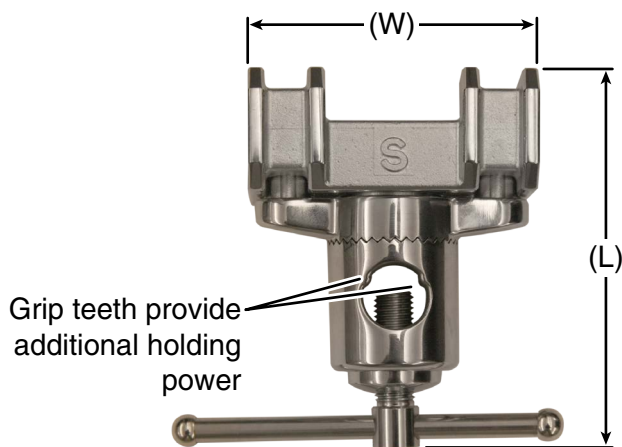

Authorized Distributors:

U.S. - SKYTRON

5085 Corporate Exchange Blvd. S.E.
Grand Rapids, MI 49512

EUROPE - SKYTRON EUROPE BV

FLORESSTRAAT 52
8022 AD ZWOLLE
THE NETHERLANDS

EC

REP



Manufactured by:

Mizuho IKA Co., Ltd.
Tokyo, Japan



TEC-D-0052 REV0
2/11